



Gynaecological chair

FG-04



The product was made of materials that possess antimicrobial/antibacterial properties.



Medisinsk utstyr og helseartikler

Telefon: 40 00 35 22
E-post: post@gymo.no
Web: www.gymo.no



Gynaecological chair

FG-04

Purpose

The FG-04 gynaecological chair has been designed to support patients during gynaecological and urological procedures.

FG-04 with foot rests WF-36.0



FG-04 with hand rest together with knee rest WF-42.0



FG-04 with leg plate WF-22.3



Standard construction

- The chair is made of varnished steel profiles covered with ABS covers
- The chair has a movable base with a central locking system
- Height adjustment, independent seat inclination and back segment inclination, simultaneous with seat segment inclination, are performed by 24 V DC electric actuators activated by a foot controller
- Trendelenburg and reverse Trendelenburg positions, simultaneous with back segment inclination and seat segment inclination, are performed by 24 V DC electric actuators activated by a foot controller
- The chair is equipped with: seamless upholstery, stainless steel gynaecological bowl, paper roll holder placed under the back segment, medical rail for additional equipment on the back rest and seat segments

Advantages of the FG-04

- Ensures patient's safety
- Easy positioning due to electric actuators
- Wide range of upholstery colours
- Seamless upholstery and ABS covers facilitate keeping the chair clean
- Applied materials are resistant to disinfecting agents
- Wide range of equipment allows the user to choose optimal version of the chair
- The chair wheels facilitate displacing the chair in the cabinet

Options

- Gynaecological chair FG-04 with foot and hand controller with programmable settings for any three seat positions

Additional equipment


- Hand rest together with knee rest WF-42.0
- Foot rest WF-36.0 (without the arm's length regulation)
- Foot stool with laminated board WF-03.5
- Foot stool with upholstered board WF-03.3
- Leg plate WF-22.3
- Foot stirrup WF-23.0
- Drip bottle holder WK-01.5
- Angle arm rest WS-07.5
- Urological bowl WF-52.2
- Wired remote controller. A programmable version is optional



Wire remote control with programmable selected items



Urological bowl WF-52.2 (option)

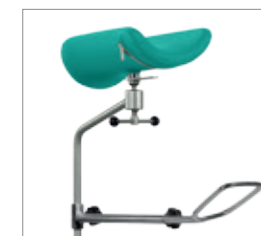

Foot stirrup WF-23.0 (option) installed on a knee rest support



Central locking system



Kitbag for paper sheet


Technical data

Length in horizontal seat and back section position	
• with leg plate	1750 mm
• without leg plate	1300 mm
Total seat and back section width	745 mm
Seat width	570 mm
Back rest width	565 mm
Seat height adjustment (in chair position, when seat segment is horizontal)	560 to 860 mm
Seat and back section height adjustment (in horizontal position)	780 to 1080 mm
Back rest inclination	-10° to +65°
Smooth seat segment inclination	-65° to +10°
Trendelenburg position	10°
Reverse Trendelenburg position	65°
Power supply	230/127/110/100 V ~, 50/60 Hz
Maximum power consumption	300 VA
Electric protection class	II
Type of application part	B
Environment influence protection degree	IP X4
Maximum permissible load	200 kg
Gynaecological bowl capacity	4l




Functions executed with use of foot controller and hand controller (optionally)


-10° to +65°

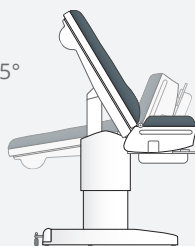


560 to 860 mm
with seat in
horizontal position

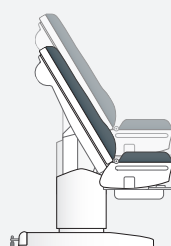

-65° to +10°

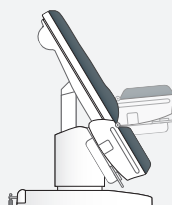

10°

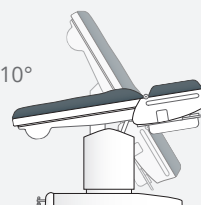

65°



The product meets the requirements of the European Directive MDD 93/42/EEC concerning safety of medical devices

Edition 05/2016/01. Famed Żywiec reserves the right to make modifications in the products in the course of technical progress. All illustrations and photos, used in this material, are for illustration and may not accurately depict the actual product.



Famed Żywiec Sp. z o.o.
ul. Fabryczna 1
34-300 Żywiec, Poland

Head Office:

tel.: +48 33 866 62 00
fax: +48 33 475 58 90

Export Department:

tel.: +48 33 866 62 58
tel.: +48 33 866 62 65
tel.: +48 33 866 62 66
fax: +48 33 861 46 78

export@famed.com.pl
www.famed.com.pl



Medisinsk utstyr og helseartikler

Telefon: 40 00 35 22
E-post: post@gymo.no
Web: www.gymo.no

