

Zen-Olite[™] with RRT

LIGHTWEIGHT PORTABLE OXYGEN CONCENTRATOR

MOBILE OXYGEN ON DEMAND FOR PATIENTS WITH ACTIVE LIFESTYLES



GIVE YOUR PATIENTS THE POWER TO STAY ACTIVE


With the Zen-O lite[™] portable oxygen concentrator from GCE Healthcare you can help your patients on long term oxygen therapy lead an independent and active lifestyle.

Weighing just 2.5kg (5.5lbs) the compact Zen-O lite™ is so portable and lightweight that there's no limit to where your patients can go.



Zen-O lite™ features a breath detection indicator

Zen-O lite™ is ultra reliable and easy to operate with a digital display, and comes with a range of additional kit that ensures patients can carry on with their lives unrestricted. Its robust and rugged design makes it suitable for everyday use.

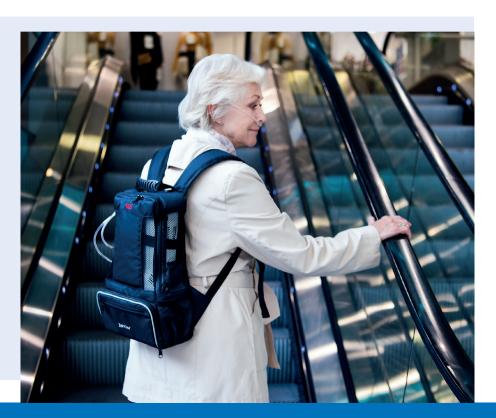

COMPREHENSIVE PACKAGE



Zen-O lite™ offers everything patients need to be independent and continue their activities outside the home.





The package includes a battery, DC and AC power supply cables, and an attractive carry bag, with an option to buy a backpack if preferred.

Also included are a user manual and quick reference guide that are easy to understand.

MEETING PATIENTS' NEEDS

Oxygen providers will benefit because Zen-O lite™ helps you to meet your patients' desire to live their lives unrestricted by respiratory disorder for as long as they are able.

The equipment's lightness, portability and flowrate adjustment technology means they can get out of their homes and carry on with the activities they love.

Rate Responsive Therapy™ (RRT)


Zen-O lite™ delivers oxygen on demand, which is triggered by the patient's inhalation.

Whether a patient is climbing the stairs or taking a walk, Zen-O lite™ will automatically respond to the breath rate by adjusting the amount of oxygen needed to ensure the patient remains saturated.

It offers an oxygen output of up to 1050ml per minute and holds one battery that can last up to 4 hours at setting 2.



BUSINESS UPLIFT FOR HOME OXYGEN **PROVIDERS**

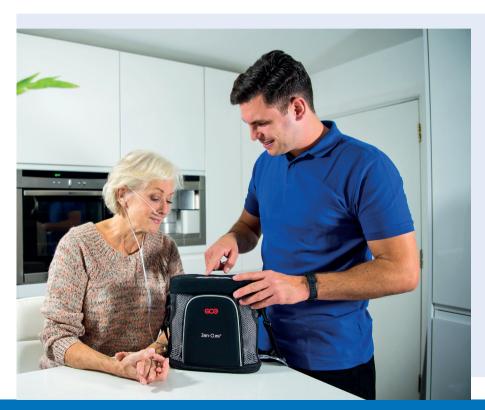

While 7en-O lite™ reduces the limitations on your patients' lifestyles, home oxygen providers will find the device is hassle-free and easy to maintain.

This ensures the cost of ownership is low, so your business success will be greater.

FIELD REPLACEABLE SIEVE MODULES

If the sieve modules need replacing, there's no need to send the device back to GCF or an authorised service centre, risking delay to the patient and the expense of supplying a new or temporary device while they wait.

It's easy to organise replacement sieve modules thanks to our robust supply chain. Then your representative or the patient simply replaces the sieve modules in a few minutes using a bespoke tool supplied.

Once the sieve module has been successfully replaced, Zen-O lite™ is easily updated from the display panel.



RELIABILITY BUILT IN

Zen-O lite™ is designed to withstand daily use and tested to independent, international standards. It is manufactured using durable components and accessories that are rigorously tested and is so quiet that your patients will sometimes forget it's there.

Zen-O lite[™] comes with a three year device warranty and sieve modules that have a oneyear warranty.

Eco Mode

The Eco Mode feature allows users to switch the operation of the device from delivery a fixed pulse of oxygen at each inhalation to a fixed volume of oxygen per minute, to allow for longer battery duration.

Auto Mode

The Auto Mode feature activates after 60. seconds if no breath is detected, the device will automatically deliver pulses at a rate of 18 breaths per minute to the user. This feature allows users to continue receiving some oxygen if the cannula is dislodged.









Zen-O lite™ has a noise level of 37dB(A) and is independently assessed as the quietest portable oxygen concentrator in its class.

YOUR BENEFITS

- > Meet the needs of active patients
- > Low cost of ownership
- > Minimal maintenance
- > Robust design fewer patient calls
- > Quick response to orders
- > Plentiful stock no waiting

YOUR PATIENTS' BENEFITS

- > More freedom to move about
- > Reliable for peace of mind
- > Lightweight, modern construction
- > Clear LCD panel
- > Easy to understand instructions
- Attractive carry case
- > Optional backpack
- > Quiet 37 dB(A)

PRODUCT CODES

ITEM NO.	DESCRIPTION
RS-00608-G-S	Zen-O lite™ portable concentrator with one 8 cell battery
RS-00608-G-D	Zen-O lite™ portable concentrator with two 8 cell batteries
RS-00601	Zen-O lite™ rechargeable battery
RS-00602	Zen-O lite™ AC power supply w/EU cord
RS-00603	Zen-O lite [™] AC power supply w/UK cord
RS-00604	Zen-O lite™ AC power supply w/US cord
RS-00605	Zen-O lite [™] DC power supply
RS-00606	Zen-O lite™ carry bag
RS-00616	Replacement sieve modules
RS-00617	Zen-O lite™ cannula wrench
RS-00619	Zen-O lite™ backpack
RS-00512	Zen-O lite™ cannula filter pk of 10
RS-00515	Zen-O™ external battery charger – US
RS-00516	Zen-O™ external battery charger – EU
RS-00517	Zen-O™ external battery charger – UK
RS-00523	Accessories bag

Each POC includes an oxygen concentrator with carry bag, battery, user manual, AC and DC power supply cords.



TECHNICAL SPECIFICATIONS

Size (W×D×H):	249 mm × 97 mm × 235 mm (9.8" × 3.8" ×9.25")
Weight:	2.5 kg (5.5 lbs) without carry bag
Power requirements:	AC adaptor: 100-240V AC (+/- 10%) 50-60 Hz in, 24V DC, 5.0A out
	DC adaptor: 11.5 - 16V DC in, 24V, 5.0A out
Purity:	87% - 96% at all settings
Maximum oxygen discharge pressure:	20.5 psi
Inspiratory trigger sensitivity:	-0.12cm/H ₂ 0
Humidity range:	5% to 93% ± 2% non-condensing
Operating altitude:	0 to 4000 meter (0' to 13000') relative to sea level, 1060 down to 575 mbar
Temperature:	
Operation:	5°C (41°F) and 40°C (104°F)
Storage:	-20°C (-4°F) and 60°C (140°F)
Setting:	Adjustable in 0.5 increments from 1.0 to 5.0
Noise level:	37 dB(A)*
Battery duration:	Approx. 4 hours*

*At setting 2.

PRODUCT CHOICE

For patients that need more oxygen, GCE Healthcare offers the Zen-O™. It's a 4,66 kg (10 lbs) portable oxygen concentrator with a higher output of two litres of oxygen per minute.

Zen-O™ is one of the few compact oxygen therapy solutions on the market to give patients the choice of either pulse or continuous flow. It means your patient can choose to have oxygen delivered constantly (continuous), such as at night when they're sleeping, as well as on demand (pulse).

QUALITY MANUFACTURING

Both Zen-O™ and Zen-O lite™ meets the requirements of the:

- > European Medical Device Directive
- > United States Food and Drug Admnistration



GCE FACILITIES



VISIT: WWW.GCEGROUP.COM

REGIONAL OFFICES

FUROPE

CZECH REPUBLIC

CZ-583 01 Chotebor Phone: +420 569 661122, Fax: +420 569 661107 sales.cz@gcegroup.com

POLAND GCE Sp z o.o.

UI. Drapińska 12 03-581 Warszawa Phone: +48 22 511 23 57, Fax: +48 22 677 70 90 sales.pl@gcegroup.com

UNITED KINGDOM & IRELAND

GCE Ltd. 100 Empress Park, Penny Lane Haydock, St Helens WA11 9DB Phone: +44 (0)1942 29 29 50, Fax: +44 (0)1942 29 29 77 sales.gb@cgeroup.com GCE S.A.S.

70 Rue du Puits Charles, BP N* 40110 FR-58403 La Charité-sur-Loire Phone: +33/3 86 69 46 00, Fax: +33/3 86 70 09 15 sales.fr@gegroup.com

PORTUGAL

GCE Portugal Rua do FeiraPark 50, Piso 2, sala 8 4520-632 São Jaão de Vêr Phone: +351 256 373 682, Fax: +351 256 378 260 sales.pt@gcegroup.com GERMANY GCE GmbH

Weyherser Weg 8 36043 Fulda Phone: +49 (0)661-8393-0, Fax: +49 661 8393 25 sales.de@cceroup.com

ROMANIA

GCE Romania S.R.L 22 Al.Puskin Street, Ap.1 Bucharest 1, 011996 Phone: +40 21 316 76 72, Fax: +40 21 316 76 74 sales.ro@qcegroup.com HUNGARY

GCE Hungária Kft.
II. Rákóczi Ferenc út 68
H-2314 Halásztelek
Phone: +36 (24) 521 200, Fax: +36 (24) 521 201
sales.hu@gcegroup.com

SPAIN

GCE Ibérica, S.L.U. Avda. de la Democracia 7–Of 311 ES-28031 Madrid Phone: +34 91 571 1470, Fax: +34 91 571 2756 sales.es@ccegroup.com TALY

GCE Mujelli spa Via F. Cervi 11 37036 San Martino B.A. (VR) Phone: +39 045 878 0 525, Fax: +39 045 878 0 750 sales.it@gcegroup.com

SWEDEN

GCE Norden AB Postal address: Box 21044, 200 21 Malmö Visitors address: Kaliwatengatan 9, 212 23 Malmö Phome: +46 (0)40 388350, Fax +46 (0)40 388353 sales.se@gcegroup.com

AMERICAS

LATIN AMERICA GCE Latín America

Po.Box: 0843-01211
Panamá, Republica de Panamá
Phone: +507 317 61 68, Fax: +507 317 65 00
sales.pa@gcegroup.com

MEXICO

GCE Gas Control Equipment SA de CV. Miguel de Cervantes Saavedra 193 Ampliación Granada Miguel, Hidalgo Distrito Federal 11529, Mexico Phone: +52 55 2626 1439 sales.pa@gcegroup.com

Gas Control Equipment Inc. 1508 Windsor Forest Trail Keller, Texas TX 76262 Phone: +1 817 268 9220, Toll Free: +1 888 659 2102 salesUSA@gcegroup.com

ASIA

CHINA

GCE Gas Control Equipment Co. Ltd. Xiao Wan Road 318, Building No.4, Fengxian District Shanghai 201401, P.R. China Phone: +86 21 37198408, Fax: +86 21 37198617 sales.cn@qcegroup.com INDIA

GCE India Pvt. Ltd.
Millers Road 59, 1st floor, Benson Town
Bangalore 560046, India
Phone: +91 802 363 1685
gce.india@gcegroup.com

RUSSIA (MOSCOW)

GCE Krass OOO Urzhumskaya Street 4 129343 Moscow, Russian Federation Phone: +7 495 745 26 99, Fax: +7 495 745 26 90 sales.ru@gcegroup.com RUSSIA (ST PETERSBURG)

GCE Krass OOO Kantemirovskaya 12A 194100 St Petersburg, Russian Federation Phone: +7 812 323 86 21, Fax: +7 812 323 86 49 sales.ru@gcegroup.com



www.gcehealthcare.com www.gcegroup.com

GCE Healthcare is part of GCE Group, a world leader in gas control equipment.

Art. Nr. DL-00498. Edition 2/2020. Alterations are subject to change without notice. 230720kk © 2020